

Iron azobenzene tetracarboxylic,



METALS • INORGANICS • ORGANOMETALLICS • CATALYSTS • LIGANDS • NANOMATERIALS • CUSTOM SYNTHESIS • cGMP FACILITIES

Catalog # 26-3725: Iron azobenzene tetracarboxylic, Porous [PCN-250(Fe)], AYRSORB™ F250 Unit sizes: 500mg, 2g and 10g

Sold in collaboration with framergy for research purposes only. US 8,940,392 B2.

Iron azobenzene tetracarboxylic (PCN-250)

Metal-organic frameworks (MOFs) are based on coordination bonds which are typically less stable than covalently bonded materials. 1-4 Therefore, the stability of MOFs is O. an important challenge to overcome. 5-15 Stronger MOFs are moisture resistant and possess a high gas storage capacity, as well as an ability to easily be scaled up. These characteristics make them more suitable for commercial sale. The problem with the manufacturing process of these MOFs is that they were made using 'one-pot' synthesis which is difficult to control¹⁶. As a way of alleviating this

challenge, an approach based on kinetic and thermodynamic control of MOF crystallization has been used to synthesize Fe-MOF single crystals with pre-synthesized metal building blocks [Fe₂ $M(\mu_3 - 0)(CH_3 COO)_6$] (M=Fe^{2+, 3+}, Co²⁺, Ni²⁺, Mn²⁺, Zn²⁺). Our PCN-250 was among the MOFs synthesize and has shown a high total H₂ volumetric uptake of 60g L-1 at 40 bar and 77 K and high total $\mathrm{CH_4}$ uptakes of 200V STP/V at 35 bar and 298 K 17,18 .

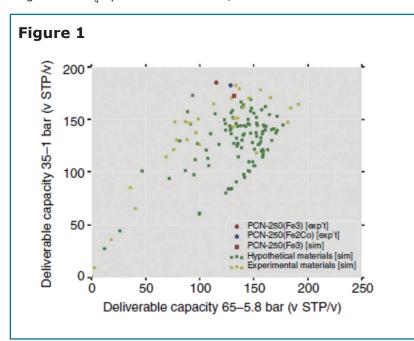


Figure 1: Performance plot of deliverable capacities using two different pressure swings. Shown are the predictions for the Fe-MOFs synthesized here (yellow) and in silico predicted structures (green). Red and blue points show PCN-250.

High volumetric gas uptake

The high methane uptake of PCN-250 is mainly due to the structure of this MOF and is a main feature of this material 18. Figure 1 below shows a performance plot of PCN-250's deliverable capacities. PCN-250 has the potential to be an excellent absorbent for methane storage if a pressure range between 1-35 bar is used with the highest methane loading at 35 bar (200 v/v and 298 K). This product also has a recorded high of H₂ uptake of 3.07 wt% and 28 g L⁻¹ at 1.2 bar and 77 K as well as a total volumetric H₂ uptake of 60 g L-1, at high pressure due to the MOFs high crystal density.

Figure 2 displays how the ligand covers each side of the cubes in PCN-250 with charged open metal sites stationed around the channels between each cube. This allows the inter space to be available for gas adsorption, which in this product has a strong interaction with both H₂ and CH₄ molecules. In addition, through inducing polarization of gas molecules through charge-induced dipole interaction a highly-efficient utilization of space for elevated volumetric gas uptake can be accomplished.

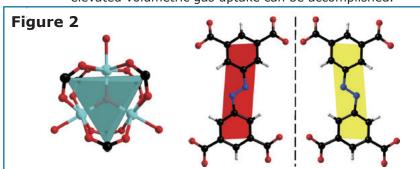


Figure 2: Structures of PCN-250

Visit www.strem.com for new product information and a searchable catalog.

Strem Chemicals, Inc. 7 Mulliken Way Newburyport, MA 01950 U.S.A

Tel: 978.499.1600 Fax: 978.465.3104 Email: info@strem.com Strem Chemicals, Inc. 15, rue de l'Atome

Zone Industrielle 67800 BISCHHEIM France Tel: (33) 03 88 62 52 60 Fax: (33) 03 88 62 26 81 Email: info.europe@strem.com

Strem Chemicals, Inc. Postfach 1215 77672 KEHL Germany Tel: 0 78 51/ 7 58 79

Email: info.europe@strem.com

Strem Chemicals UK Ltd.
An Independent Distributor of Strem Chemicals Products Newton Hall, Town Street Newton, Cambridge England CB22 77F Tel: +44 (0)1223 873 028

Fax: +44 (0)1223 870207 Email: enquiries@strem.co.uk

0718

Figure 3

Figure 3: Structure drawing of PCN-250

Chemical stability of PCN-250

Studies have shown that the chemical stability of PCN-250 remains even if one iron atom is substituted for a softer Lewis acid metal, M(II), in the μ_3 -oxo cluster. The product remains intact in a range of pH values scaling from 1 to 11 over a period of 24 hours. It was also found that PCN-250 remained strong in H_2 O after 6 months in neutral conditions.

After multiple pH treatments, there was barely any change in the material's N_2 adsorption isotherms leading researchers to conclude that no MOF decomposition or phase transition occurred.

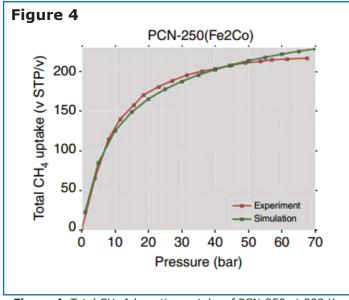


Figure 4: Total CH₄ Adsorption uptake of PCN-250 at 298 K.

PCN-250's high gas uptake and exceptional stability and scalability make it an auspicious material for natural gas storage for a wide range of applications including power systems and gas purification.

Key Properties

- · Color and form: dark red-brown pwdr.
- · Superior uptake of hydrogen and methane
- · Stable in water and aqueous solutions

References

- 1. Zhou, H.-C., Yaghi, O. M. & Long, J. R. Introduction to metal-organic frameworks. Chem. Rev. 112, 673-674 (2012).
- 2. Fe'rey, G. et al. A chromium terephthalate-based solid with unusually large pore volumes and surface area. Science 309, 2040–2042 (2005).
- 3. Serre, C. et al. Role of solvent-host interactions that lead to very large swelling of hybrid frameworks. Science 315, 1828-1831 (2007).
- 4. Banerjee, R. et al. High-throughput synthesis of zeolitic imidazolate frameworks and application to CO₂ capture. Science 319, 939–943 (2008).
- 5. Kitagawa, S., Kitaura, R. & Noro, S. Functional porous coordination polymers. Angew. Chem. Int. Ed. 43, 2334–2375 (2004).
- 6. Herm, Z. R. et al. Separation of hexane isomers in a metal-organic framework with triangular channels. *Science* **340**, 960–964 (2013).
- 7. Deng, H. et al. Large-pore apertures in a series of metal-organic frameworks. Science **336**, 1021–1023 (2012).
- 8. Bloch, E. D. et al. Hydrocarbon separations in a metal-organic framework with open iron (II) coordination sites. Science 335, 1606–1610 (2012).
- 9. Farha, O. K. et al. De novo synthesis of a metal-organic framework material featuring ultrahigh surface area and gas storage capacities. *Nat. Chem.* 2, 944–948 (2010).
- Cavka, J. H. et al. A new zirconium inorganic building brick forming metal organic frameworks with exceptional stability. J. Am. Chem. Soc. 130, 13850–13851 (2008).
- 11. Fe'rey, G. et al. A hybrid solid with giant pores prepared by a combination of targeted chemistry, simulation, and powder diffraction. *Angew. Chem. Int. Ed.* **43**, 6296–6301 (2004).
- 12. Surble', S. et al. A new isoreticular class of metal-organic-frameworks with the MIL-88 topology. Chem. Commun. 3, 284–286 (2006).
- 13. Dincă, M. et al. Hydrogen storage in a microporous metal—organic framework with exposed Mn2t coordination sites. *J. Am. Chem. Soc.* **128**, 16876–16883 (2006).
- 14. Li, J.-R., Sculley, J. & Zhou, H.-C. Metal-organic frameworks for separations. Chem. Rev. 112, 869-932 (2012).
- 15. Yoon, M., Srirambalaji, R. & Kim, K. Homochiral metal-organic frameworks for asymmetric heterogeneous catalysis. Chem. Rev. 112, 1196–1231 (2012).
- 16. Tranchemontagne, D. J., Mendoza-Corte's, J. L., O'Keeffe, M. & Yaghi, O. M. Secondary building units, nets and bonding in the chemistry of metal—organic frameworks. *Chem. Soc. Rev.* 38, 1257–1283 (2009).
- 17. Murray, L. J., Dincă, M. & Long, J. R. Hydrogen storage in metal-organic frameworks. Chem. Soc. Rev. 38, 1294-1314 (2009).
- Makal, T. A., Li, J.-R., Lu, W. & Zhou, H.-C. Methane storage in advanced porous materials. Chem. Soc. Rev. 41, 7761–7779 (2012).
- Schoedel, A. & Zaworotko, M. J. [M3(μ3-O)(O2CR)6] and related trigonal prisms: versatile molecular building blocks for crystal engineering of metal-organic material platforms. Chem. Sci. 5, 1269–1282 (2014).
- 20. Liu, Y. et al. Assembly of metal—organic frameworks (MOFs) based on indiumtrimer building blocks: a porous MOF with soc topology and high hydrogen storage. *Angew. Chem. Int. Ed.* **46**, 3278–3283 (2007).

Please see our Chemiker XXIX (2016) for a detailed article about this product.

Visit www.strem.com for new product information and a searchable catalog

Strem Chemicals, Inc. 7 Mulliken Way Newburyport, MA 01950 U.S.A Tel: 978.499.1600

Tel: 978.499.1600 Fax: 978.465.3104 Email: info@strem.com Strem Chemicals, Inc. 15, rue de l'Atome Zone Industrielle 67800 BISCHHEIM France

Tel: (33) 03 88 62 52 60 Fax: (33) 03 88 62 26 81 Email: info.europe@strem.com Strem Chemicals, Inc. Postfach 1215 77672 KEHL Germany Tel: 0 78 51/ 7 58 79

Email: info.europe@strem.com

Strem Chemicals UK Ltd. An Independent Distributor of Strem Chemicals Products Newton Hall, Town Street

Newton Hall, Town Street Newton, Cambridge England CB22 7ZE Tel: +44 (0)1223 873 028

Fax: +44 (0)1223 870207 Email: enquiries@strem.co.uk